



Field Service Procedure

Part Number: SP00277

Rev: A

Date: 18 November 2003

© 2003 Draeger Medical, Inc.

Flip Up Tray Installation on Fabius Tiro™

[RETURN TO SERVICE PROCEDURE TABLE OF CONTENTS](#)
[RETURN TO CD-ROM TABLE OF CONTENTS](#)

Flip-up Tray Mounting Instructions

Purpose: The purpose of this document is to provide instruction for the installation and operation of the Flip-up Tray.

Parts Reference

The following parts and hardware are included with this installation kit.

Item #	Description	Qty
1	Flip up tray	1
2	GCX handle	1
3	Resuscitator mounting hook	1
4	#10 self tapping screw	2
5	5/32" Tension wrench	1

INSTALLATION PROCEDURE:

1. Remove the two M8 socket head cap screws (item 1), two spacers (item 2), and Draeger rail (item 3) from the side of the trolley (Figure 1).

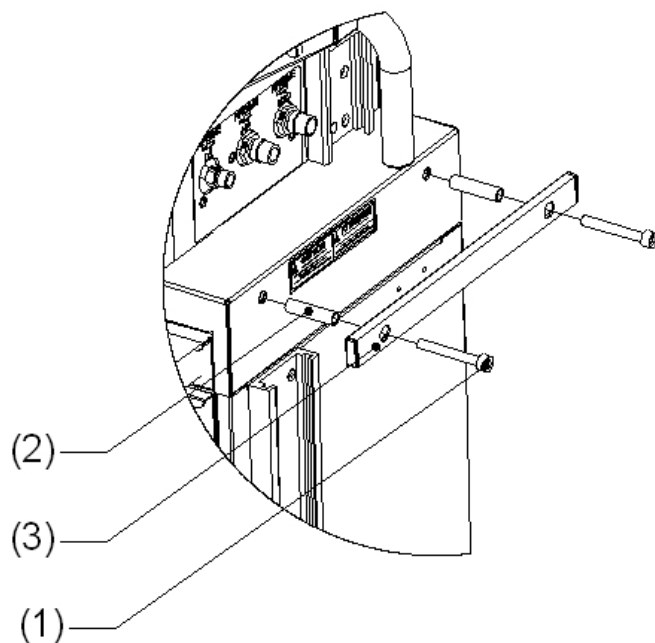


Figure 1. Remove Draeger Rail

SP00277
10/6/03

2. Insert the Flip-up Tray slide (item 1) in the GCX channel (Figure 2).
3. Move the Flip-up Tray to the desired height and tighten the two clamping screws (item 2) on slide bottom to secure the Flip-up Tray in the channel (Figure 2).

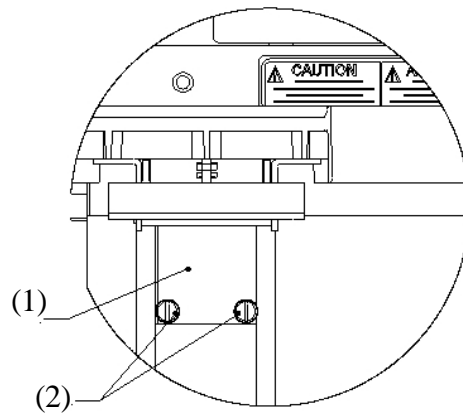


Figure 2. Attach Flip-up Tray

4. Insert the handle (item 1) in the GCX channel (item 2) as shown in Figure 3.
5. Route the hoses (item 3) under the handle, then adjust handle until handle rests on top surface of trolley (Figure 3).
6. Tighten set screws (item 4) on the handle with wrench provided.

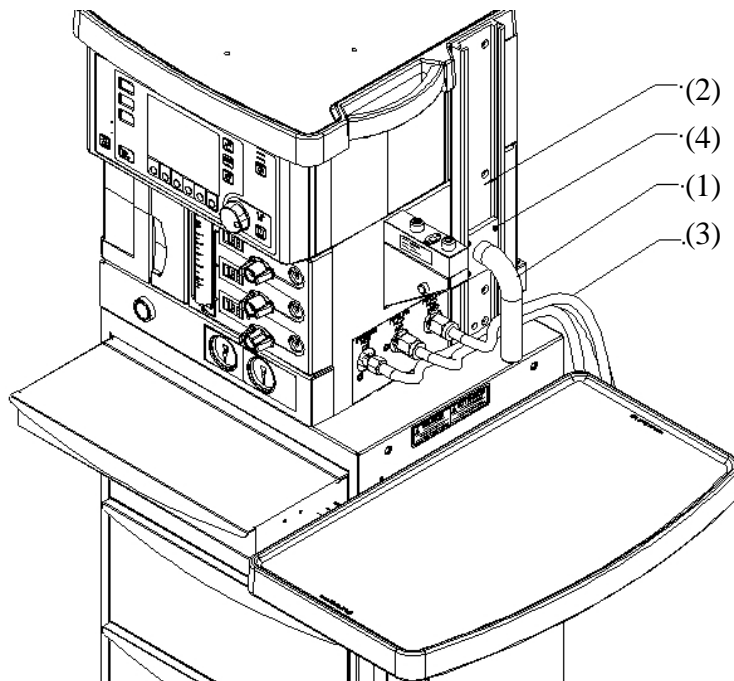


Figure 3. Hose routing and handle attachment

7. Attach resuscitator hook (item 1) to back of Fabius Tiro using two #10 self tapping screws (item 2) provided (Figure 4).

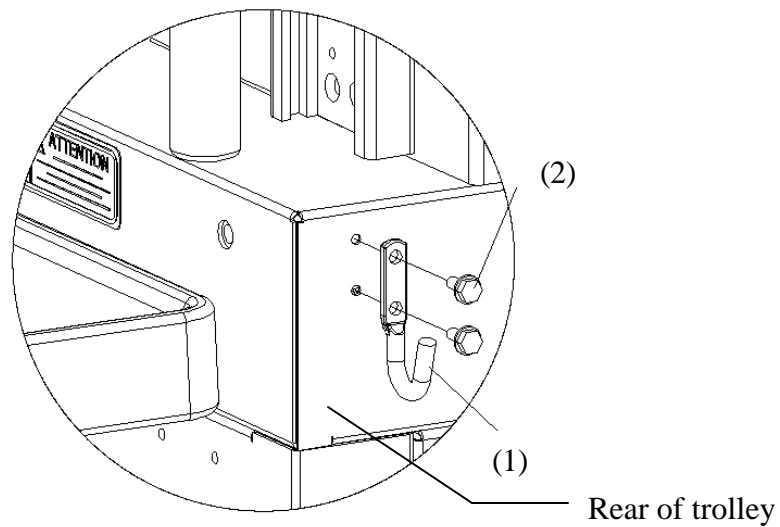


Figure 4. Attach resuscitator hook

SP00277
10/6/03

OPERATION OF FLIP UP TRAY:

1. To open tray, lift tray until it locks in position (locking is identified when click is heard).
2. To release tray, pull the tray handle (item 1) forward (Figure 5).

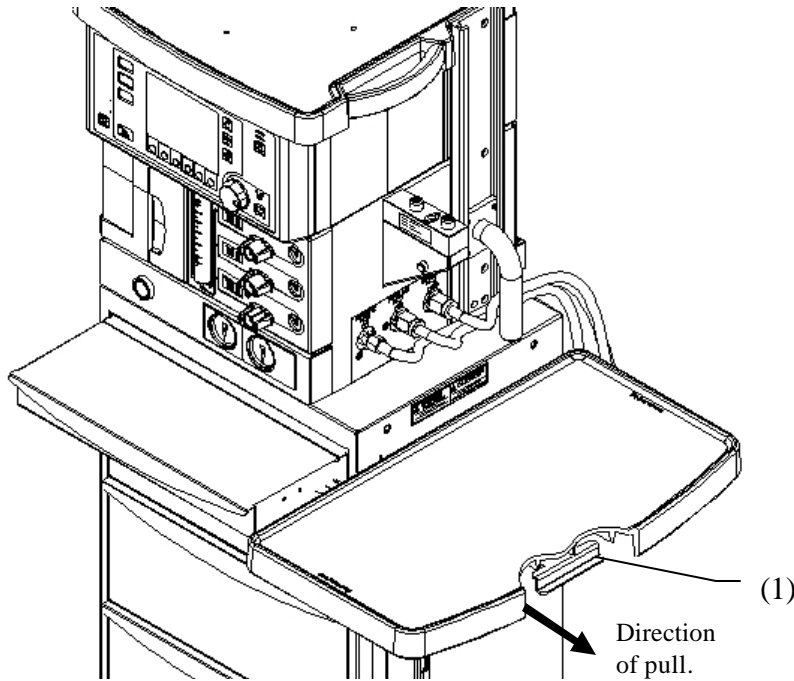


Figure 5. Tray in extended position

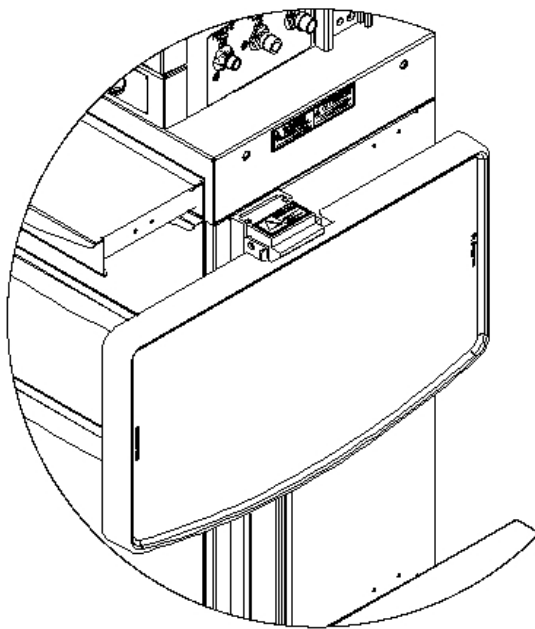


Figure 6. Tray in folded position

[RETURN TO SERVICE PROCEDURE TABLE OF CONTENTS](#)
[RETURN TO CD-ROM TABLE OF CONTENTS](#)



DrägerService is a division of
Draeger Medical, Inc.
3122 Commerce Drive
Telford, PA 18969
Tel: (215) 721-5402
(800) 543-5047
Fax: (215) 721-5784
Web: www.draegermedical.com
Printed in the U.S.A.